



Product	Description
---------	-------------



FCSLWAF21	21mm, dual safety, folding anvil, 4.5mm staple height Extra Long
FCSLWAF25	25mm, dual safety, folding anvil, 4.8mm staple height, Extra Long
FCSLWAF29	29mm, dual safety, folding anvil, 4.8mm staple height, Extra Long
FCSLWAF32	32mm, dual safety, folding anvil, 5.0mm staple height, Extra Long




FCSME21	21mm, dual safety, folding anvil, 4.5mm staple height
FCSME24	24mm, dual safety, folding anvil, 4.5mm staple height
FCSME26	26mm, dual safety, folding anvil, 4.8mm staple height
FCSME29	29mm, dual safety, folding anvil, 4.8mm staple height
FCSME32	32mm, dual safety, folding anvil, 5.0mm staple height

Laparoscopic circular staplers

A FULL RANGE OF ADVANCED CIRCULAR STAPLERS TO COVER ALL SURGEONS NEEDS




Standard
&
Extra long

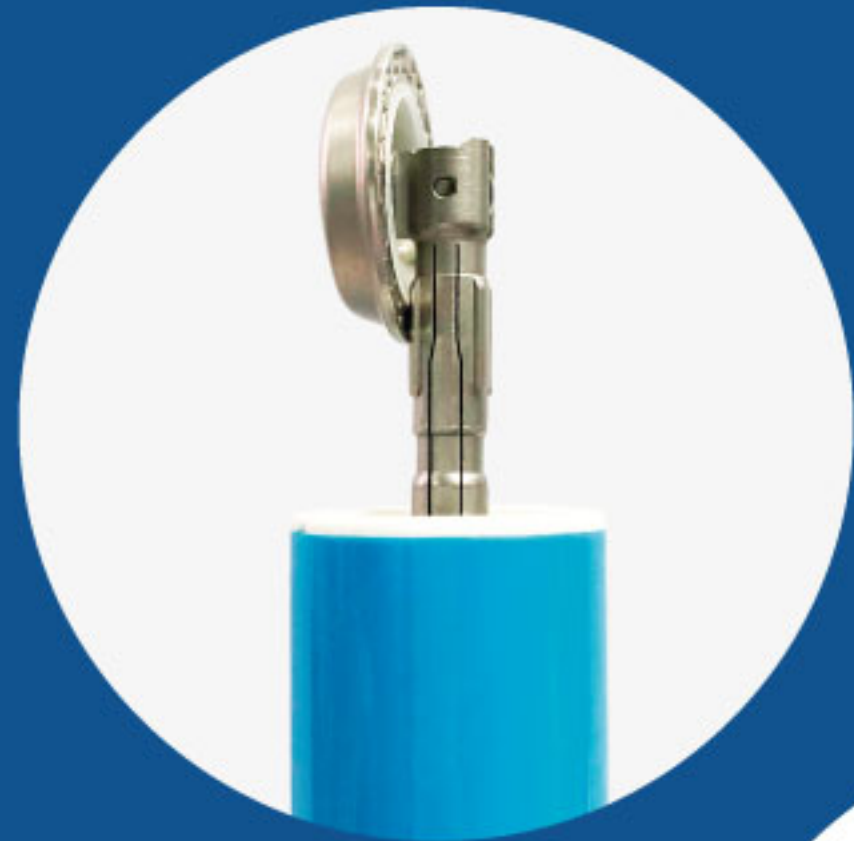

Unique 90 degrees Folding anvil
For easy & safe removal through anastomosis,
minimizing potential tissue trauma




Improved cutting washer
Gives surgeon a clear feedback of firing

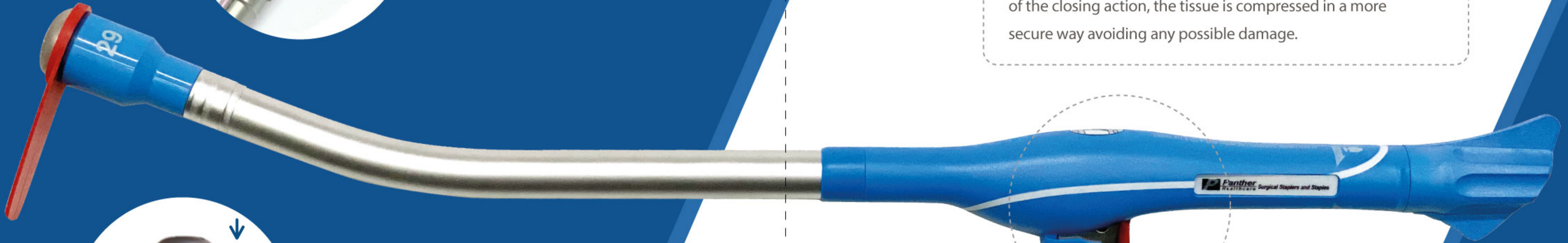

Cartridge and Anvil Self-Lock
Prevents slipping out of the tissue and
eliminates the problem when the anvil
returns back during the firing action.

Hermetic design
Hermetic design
for laparoscopic surgery,
prevents the loss
of pneumoperitoneum

LAPAROSCOPIC CIRCULAR STAPLERS

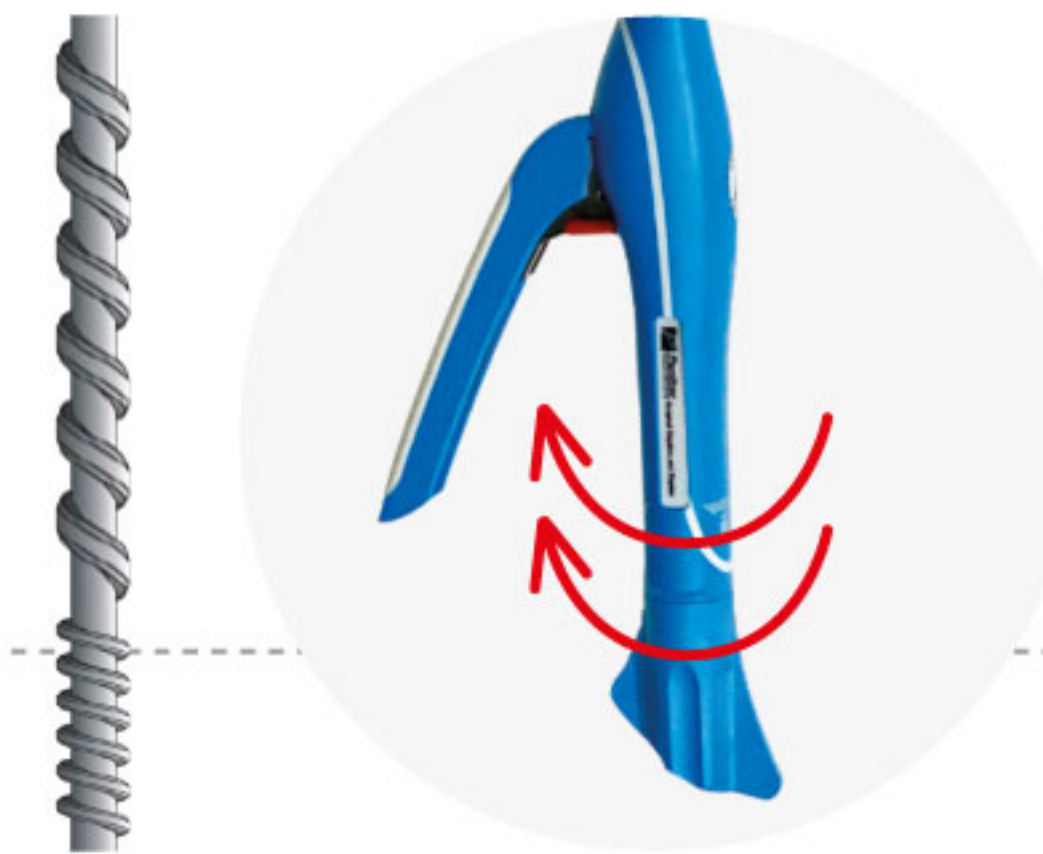

Dual-Speed Function

When closing the instrument, at the final stages of closing, the closure obtained per each turn of the closing knob is smaller than at the beginning. Provides better control of the closing action, the tissue is compressed in a more secure way avoiding any possible damage.



**Dual safety —
Patent design**
Prevents from
accidental misfiring