

Flexella Plus

SX-ELLA Stent Esophageal



Flexella Plus



Durable silicone covering

Removable esophageal stent

Anti-migration and atraumatic flares

Excellent flexibility



Basic description

SX-ELLA Stent Esophageal – Flexella Plus

Flexella Plus is designed for relieving malignant or benign esophageal strictures if standard operation is contraindicated. It can be used for esophagorespiratory fistulas treatment as well.







Stent system for fluoroscopic and endoscopic insertion.
Self-expandable fully covered nitinol stent preloaded in delivery system.

Features / Benefits

- ☑ 3 types of delivery system: PULL 16 F / 22 F, PULL 18 F / 28 F and PUSH 14 F / 28 F.
- ☑ PUSH delivery system - strictures located in upper part of esophagus.
- ☑ PULL (28 F / 18 F) delivery system with splittable olive tip - automatic detachment while stent starts releasing; olive entrapment prevention; delivery system safe removal.
- ☑ Retrieval loop made of medical grade alloy available in two options: on the proximal end for stent diameter 20 mm and on both stent ends for stent diameter 18 mm - high mechanical durability and acid-resistance; stent flare centripetal constriction; suitable for extraction by ELLA Extractor.
- ☑ Radiopaque markers - accurate stent positioning.
- ☑ Anti-migration design - flared design; reduced migration rates.
- ☑ Atraumatic stent ends - esophageal wall irritation and hyperplasia reduction.
- ☑ Excellent flexibility - smooth alignment and adequate stent expansion along anatomical curves.
- ☑ Durable silicone stent covering - tumor in-growth prevention, esophagorespiratory fistulas occlusion.
- ☑ Optional anti-reflux valve - gastroesophageal reflux and aspiration prevention.

Additional description

For stent diameter Ø 20 mm - PUSH delivery system (06S model) 4.7 mm / 9.4 mm (14 F / 28 F), available in standard length 55 cm.
For stent diameter Ø 20 mm - PULL delivery system (07S model) 6.0 mm / 9.4 mm (18 F / 28 F), available in standard length of 75 cm.
For stent diameter Ø 18 mm - PULL delivery system (07S model) 5.2 mm / 7.2 mm (16 F / 22 F), available in standard length of 75 cm.
Recommended guide wire: 0.035" (0.89 mm) - ultra stiff - minimum length of 220 cm.

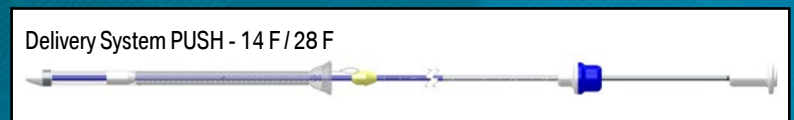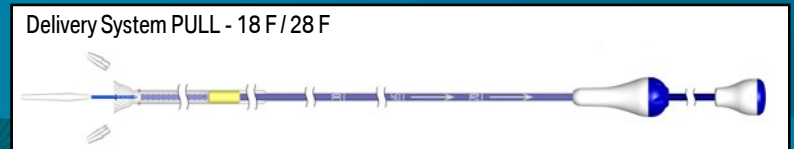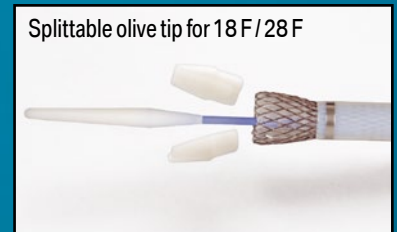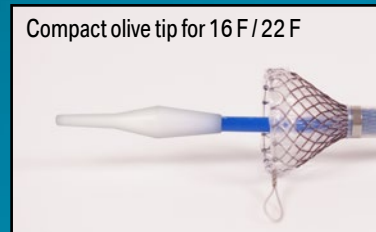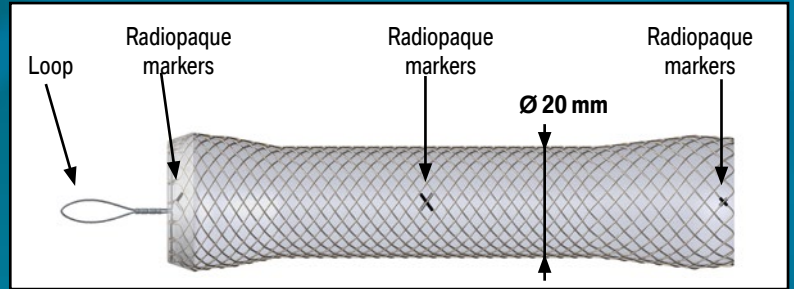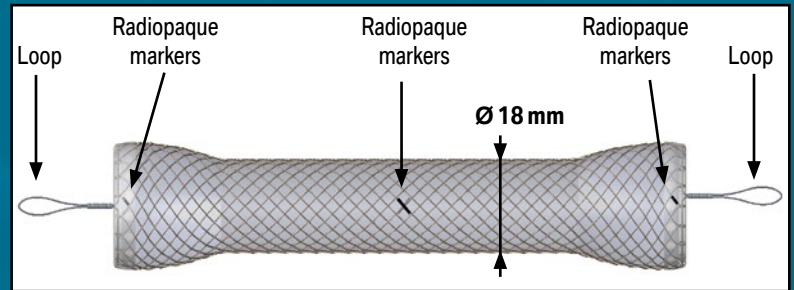

Available sizes

| REF No. | Stent diameter [mm] Proximal flare / Body / Distal flare | Stent nominal length [mm] | Anti-reflux valve | Delivery system |
|------------------|---|---------------------------|-------------------|--------------------------------|
| 019-06S-20-085 | 25/20/24 | 85 | No | 14 F / 28 F - proximal release |
| 019-06S-20-110 | 25/20/24 | 110 | No | 14 F / 28 F - proximal release |
| 019-07S-20-085 | 25/20/24 | 85 | No | 18 F / 28 F - distal release |
| 019-07S-20-110 | 25/20/24 | 110 | No | 18 F / 28 F - distal release |
| 019-07S-20-135 | 25/20/24 | 135 | No | 18 F / 28 F - distal release |
| 019-07S-20-150 | 25/20/24 | 150 | No | 18 F / 28 F - distal release |
| 019-07S-20-085-0 | 25/20/24 | 85 | Yes | 18 F / 28 F - distal release |
| 019-07S-20-110-0 | 25/20/24 | 110 | Yes | 18 F / 28 F - distal release |
| 019-07S-20-135-0 | 25/20/24 | 135 | Yes | 18 F / 28 F - distal release |
| 019-07S-20-150-0 | 25/20/24 | 150 | Yes | 18 F / 28 F - distal release |
| 019-07S-18-085 | 25/18/25 | 85 | No | 16 F / 22 F - distal release |
| 019-07S-18-110 | 25/18/25 | 110 | No | 16 F / 22 F - distal release |
| 019-07S-18-135 | 25/18/25 | 135 | No | 16 F / 22 F - distal release |
| 019-07S-18-150 | 25/18/25 | 150 | No | 16 F / 22 F - distal release |
| 019-07S-18-085-0 | 25/18/25 | 85 | Yes | 16 F / 22 F - distal release |
| 019-07S-18-110-0 | 25/18/25 | 110 | Yes | 16 F / 22 F - distal release |
| 019-07S-18-135-0 | 25/18/25 | 135 | Yes | 16 F / 22 F - distal release |
| 019-07S-18-150-0 | 25/18/25 | 150 | Yes | 16 F / 22 F - distal release |



Distributor:

Manufacturer:

ELLA-CS, s.r.o.

Milady Horakove 504/45, Trebes,
500 06 Hradec Kralove, Czech Republic

Tel.: +420 495 279 111 Fax: +420 495 265 655
E-mail: info@ellacs.eu URL: www.ellacs.eu